



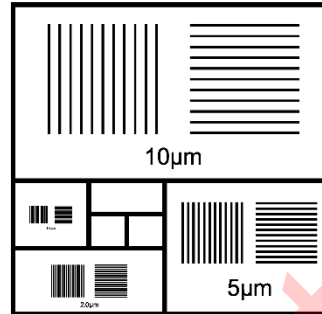

# AISthesis Products, Inc.

Advanced Imaging Products for Nanotechnology, Engineering and Life Sciences

## Certificate of Calibration for Pelcotec™ XY-Critical Dimension Standard

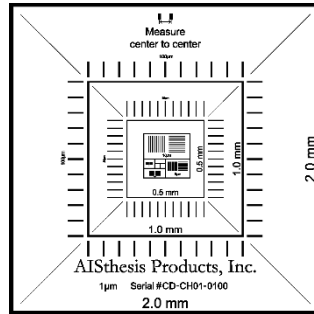

**Product Number:** Pelcotec™ 705-1 CDMS-XY-1C-Etched

**Manufactured for and Distributed by:**

**Product Description:** 2.5x2.5mm Pelcotec™ 2mm-1μm XY-Critical Dimension Magnification Standard



**Product Serial Number:** CD-CH01-1234

The accuracy of these products was determined by reference comparison to working standards traceable to the National Institute of Standards and Technology (NIST).

Line	Direction	Certified average pitch	Number of lines averaged	Certified distance between first and last line (1σ)	Total expanded uncertainty (3σ)
2.0mm	X	2.00 mm	2	2.00 mm ± 0.03%	± 0.09%
1.0mm	X	1.00 mm	2	1.00 mm ± 0.03%	± 0.09%
0.5mm	X	0.500 mm	2	0.500 mm ± 0.03%	± 0.09%
0.10mm	X	100.04 μm	11	1.00 mm ± 0.03%	± 0.09%
50μm	X	50.00 μm	11	0.500 mm ± 0.03%	± 0.09%
10μm	X	10.00 μm	11	100.04 μm ± 0.03%	± 0.09%
5μm	X	5.01 μm	12	55.10 μm ± 0.03%	± 0.09%
2μm	X	2.01 μm	16	30.13 μm ± 0.03%	± 0.09%
1μm	X	1.00 μm	17	16.07 μm ± 0.03%	± 0.09%
2.0mm	Y	2.00 mm	2	2.00 mm ± 0.03%	± 0.09%
1.0mm	Y	1.00 mm	2	1.00 mm ± 0.03%	± 0.09%
0.5mm	Y	0.500 mm	2	0.500 mm ± 0.03%	± 0.09%
0.10mm	Y	100.13 μm	11	1.00 mm ± 0.03%	± 0.09%
50μm	Y	50.00 μm	11	0.500 mm ± 0.03%	± 0.09%
10μm	Y	10.01 μm	11	100.13 μm ± 0.03%	± 0.09%
5μm	Y	5.01 μm	12	55.12 μm ± 0.03%	± 0.09%
2μm	Y	2.01 μm	16	30.15 μm ± 0.03%	± 0.09%
1μm	Y	1.01 μm	17	16.08 μm ± 0.03%	± 0.09%

The certified average pitch is derived from the stated certified length that was determined using 10 measurements (taken center-to-center) over the stated number of lines (i.e., length divided by the number of lines minus one). The total expanded uncertainty includes both Type A and Type B uncertainties corrected for sample size using an appropriate Student t-factor.

Equipment used:

Instrument	Model number	Serial #	NIST Certified CD	Resolution	Repeatability
FE-SEM	FEI Apreo 2	9958357	CD-PG01-0211	0.9nm	0.03%

D.S. Finch  
Certified by

\_\_\_\_\_  
Signature

01/26/2024  
Date

This certificate shall not be reproduced without the permission of AIStthesis Products, Inc.  
P.O. Box 1950, Clyde, North Carolina 28721 Tel: 828.627.6555 E-mail: CD@aistthesisproducts.com

Example Certificate Only