


AISThesis Products, Inc.

Advanced Imaging Products for Nanotechnology, Engineering and Life Sciences

Wafer Level Certificate of Traceability for Pelcotec™ Critical Dimension Standard

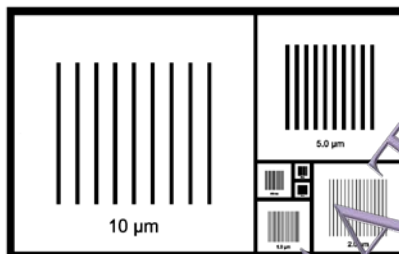

Product Number: Pelcotec™ 683-01 CDMS-0.1T

Manufactured for and distributed by:

Product Description: 2.5x2.5mm Pelcotec™ 2mm-100nm Critical Dimension Magnification Standard



Wafer Identifier: CD-XE01

The accuracy of these products was determined by reference comparison to working standards traceable to the National Institute of Standards and Technology (NIST), Test No. 861/280822-11.

Line	Average pitch of wafer	Number of lines averaged	Average pitch uniformity (1σ uncertainty)	Total expanded uncertainty (3σ) average pitch for wafer*
2.0mm	2.00 mm	2	± 2µm (±0.10%)	± 7µm (±0.35%)
1.0mm	1.00 mm	2	± 1µm (±0.10%)	± 3.5µm (±0.35%)
0.5mm	0.500 mm	2	± 0.5µm (±0.10%)	± 1.75µm (±0.35%)
0.25mm	0.250 mm	2	± 0.25µm (±0.10%)	± 0.9µm (±0.35%)
10µm	10.00 µm	9	± 0.01µm (±0.10%)	± 0.035µm (±0.35%)
5µm	5.00 µm	12	± 0.01µm (±0.20%)	± 0.035µm (±0.70%)
2µm	2.00 µm	16	± 0.004µm (±0.20%)	± 0.014µm (±0.70%)
1µm	1.00 µm	17	± 0.002µm (±0.20%)	± 0.007µm (±0.70%)
500nm	500.8 nm	20	± 1.00nm (±0.20%)	± 3.5nm (±0.70%)
250nm	250.7 nm	21	± 0.50nm (±0.20%)	± 1.75nm (±0.70%)
100nm	100.3 nm	52	± 0.20nm (±0.20%)	± 0.75nm (±0.70%)

* The 3σ uncertainty (95% confidence interval) average pitch is determined using a minimum of nine die per production wafer. Each average pitch is determined using 100+ measurements on each die averaged over the stated number of lines. The total expanded uncertainty includes both Type A and Type B uncertainties corrected for sample size using an appropriate Student t-factor.

Equipment used:

Instrument	Model number	Serial #	NIST Certified CD/Recalibration	Resolution	Repeatability
FE-SEM	FEI Verios	9922551	CD-PG01-0211/June 2018	0.9nm	0.03%

Certified by _____

Signature _____

Date _____

This certificate shall not be reproduced without the permission of AISThesis Products, Inc.
P.O. Box 1950, Clyde, North Carolina 28721 Tel: 828.627.6555 E-mail: CD@aisthesisproducts.com