


## HIGH ORDERED PYROLYTIC GRAPHITE ZYB QUALITY

Product No. 626-10


### DESCRIPTION:

Graphite consists of identical stacked planes (*Figure 1*) of covalently bonded graphene. The layers are weakly held together and can be easily exfoliated. High Ordered Pyrolytic Graphite (HOPG) is a crystallographically textured graphite with the majority of the graphene sheets stacked parallel to the surface of the material. The surface has outstanding regularity and smoothness at the nanoscale (*Figure 2*). HOPG is ideal as a test substrate for STM and AFM applications. An atomically fresh surface can be easily exposed via exfoliation with tape.

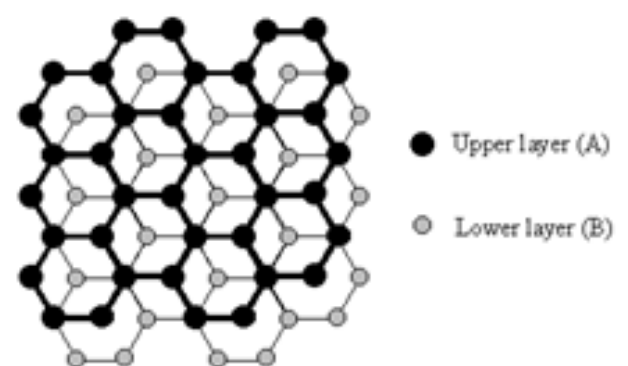

*Figure 1:*

Graphite structure can be described as an alternate succession of these plans, such as:  
“... ABABAB ...”

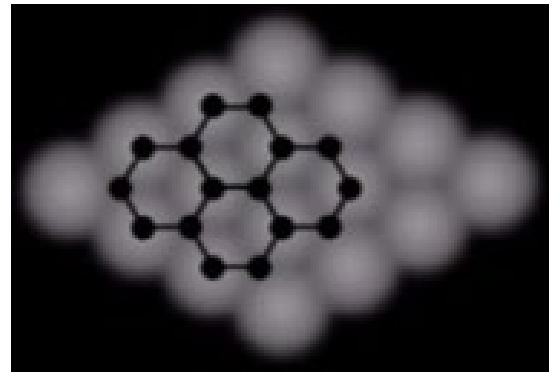

*Figure 2:*

Typical STM image of HOPG surface. Corresponding fragment of graphene structure is superimposed.

### APPLICATIONS:

- Scanning Tunneling Microscopy Calibration Samples
- SPM Substrates
- Monochromators
- X-Ray Diffractometry
- Neutron Scattering Experiments and Diffractometry

### SPECIFICATIONS:

- ZYB Quality / Mosaic Spread  $0.8^{\circ} \pm 0.2^{\circ}$
- Substrate Size 1 x 10 x 10mm
- Grain Size up to  $1\mu\text{m}$  – Double-sided

This HOPG is double-sided and both top and bottom surfaces may be used for imaging.

**TED PELLA, INC.**

*Microscopy Products for Science and Industry*

P.O. Box 492477, Redding, CA 96049-2477, U.S.A.  
Telephone: 530-243-2200; 800-237-3526 (U.S.A. or Canada) • FAX: 530-243-3761  
Email: sales@tedpella.com • Web Site: <http://www.tedpella.com>

626-10 TN V5 04052021

Page 1 of 1